



# TissueScope™ LE Slide Scanner

The best in whole slide imaging, only better.

The award-winning<sup>1</sup> TissueScope™ LE scanner sets the standard for versatility, throughput, and ease of use. Scan up to twelve 1" x 3" brightfield slides, three 2" x 3" slides, or any size slide up to 6" x 8".

## Key Features

- Easy to use
- Fast: < 30 seconds per slide<sup>2</sup>
- Sharp, clear image quality up to 40X magnification
- Scan up to twelve standard 1" x 3" slides, or any size up to 6" x 8" for unmatched versatility
- Z-Stack scanning
- Image-based barcode reading
- Non-proprietary file format, DICOM compatible
- Optimize workflow with TissueSnap™ preview station
- Images immediately available for viewing after scanning

<sup>1</sup> Validated Versatility Award, 3rd International Scanner Competition



## Technical Specifications

<b>Imaging Mode</b>	Brightfield
<b>Slide Formats</b>	25mm x 75mm (1" x 3") 50mm x 75mm (2" x 3") Any slide up to 150mm x 200mm (6" x 8")
<b>Optical Resolution</b>	0.75 NA
<b>Resolution (µm/pixel)</b>	0.25 at 40X 0.5 at 20X
<b>Scan Speed</b>	< 30 seconds per slide <sup>2</sup>
<b>File Format</b>	Non-proprietary 24-bit RGB Pyramidal BigTIFF Uncompressed, JPEG or JPEG2000 Compression Export to JPEG, Flat TIFF and LZW Compressed TIFF
<b>Dimensions &amp; Weight</b>	61cm (24.1") Width 56cm (22.2") Length 40cm (15.9") Height 39Kg (85lbs) Weight
<b>Certification</b>	CSA/UL
<b>Warranty</b>	One year

<sup>2</sup> 15mm x 15mm @ 40X magnification

## Different Slide Sizes? Difficult Slides? No Problem.


Rely on our unique tray-based slide holder to safely scan multiple slide sizes from 1" x 3" to 6" x 8".

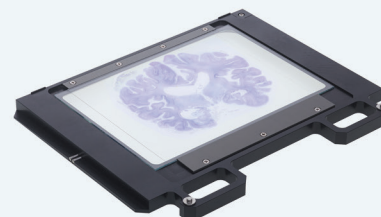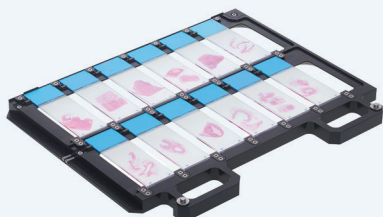

Have faint tissue, chipped slides, thick slides, or multiple bar codes? We can image them.

## Optimize your workflow with the TissueSnap™ preview station

The TissueSnap™ preview station works with our TissueScope™ LE scanners to offload time consuming tasks such as preview scanning and batch pre-processing so the scanner can work at peak throughput and performance.

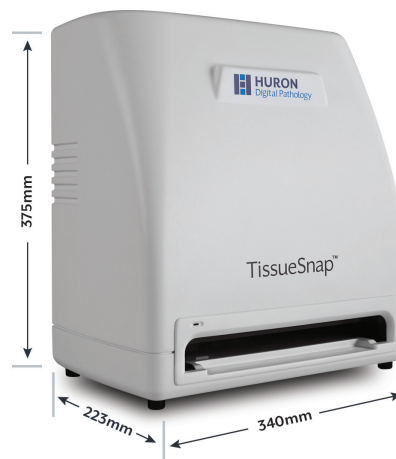

## Bundle image management and analysis



### Image Management

Connect your TissueScope™ LE scanner to Huron's SlideVault™ image management system for a complete slide scanning and management solution.

### Image Analysis

We partner with the leading image analysis software vendors to provide image analysis solutions for research. Our Lagotto™ Image Search Engine for Pathology enables the indexing and searching of whole slide images.

## About Huron Digital Pathology

We believe that innovation is the key to making digital pathology a ubiquitous reality. But for us, innovation is about more than just delivering great image quality and fast scanning speeds. It's also about designing products that are simple to use, easy to integrate with your workflow, and have attractive price-performance characteristics.

Based in St. Jacobs, Ontario, Canada, Huron Digital Pathology has a 20 year history designing sophisticated imaging instrumentation. Our end-to-end digital whole slide scanning solutions for digital pathology incorporate our award-winning TissueScope™ digital slide scanners, image management software, and our workflow-enhancing accessories, which include our innovative TissueSnap™ preview scanning station.



Manufactured by **Huron Technologies International Inc.**, of which Huron Digital Pathology is an operating division.  
1620 King St. North, PO Box 504, St. Jacobs, ON N0B 2N0 Canada

TissueScope™, TissueSnap™, Lagotto™, SlideVault™ and Huron Digital Pathology™ are trademarks of Huron Technologies International Inc.

The TissueScope™ scanner is based on one or more of the following: U.S. Patents 8,896,918 B2/8,655,043/7,218,446 B2 / 6,072,624; EU patent 2758825; U.S. and International patents pending.

Document #SMB-0004-03